SOLO^{RX}

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Table of contents

Safety regulations	3
Introduction	10
Technical specifications	11
Product description	13
Settings and adjustments	18
Usage	21
Maintenance	29
Transport and storage	33
Malfunctions	35
Warranty	36
Dispose of product and packaging	37

Information for people with visual impairments can be obtained from Life & Mobility. Contact details can be found on the last page.

Safety regulations

As a mobility scooter user, you are always responsible for complying with local safety instructions and guidelines. It's important to first gain riding experience under the supervision of an experienced instructor. Your dealer can advise you on this.

Life & Mobility accepts no liability for any loss or injury caused by negligence or failure to follow safety instructions.

Life & Mobility and the competent authority must always be informed in the event of a serious incident.

NB

- *This symbol in the manual indicates additional instructions. Always follow these instructions. Failure to do so can adversely affect the scooter's operation and may cause bodily injury or damage to the scooter or the surrounding area.*

General

Warning

- *The construction of the scooter may not be changed under any circumstances.*
- *Keep the scooter away from open flames.*
- *Before installing, maintaining, or operating any equipment, consult this manual.*

Be careful:

- *Use the basket only for transporting objects.*
- *Do not use the scooter to pull vehicles or carts, for example.*
- *Be aware that there are fixed and loose parts protruding from the scooter and its surroundings. These can cause damage to the scooter and the surrounding area.*

NB

- This scooter can transport one person.
- Use only original parts and accessories from Life & Mobility.
- The scooter features fully programmable electronic controls. Programming may only be modified by Life & Mobility.
- The safety instructions and stickers must never be covered or removed.
- Always turn off the scooter and remove the ignition key when leaving it somewhere. This prevents unauthorized use and/or theft.

User Requirements

Warning

- You may not use the scooter if you are under the influence of substances that could affect your driving ability.
- You must be physically and mentally able to operate the scooter in the intended environment.

NB

- Your eyesight must be sufficient to operate the scooter safely in any situation that may arise during its use.
- Please observe applicable traffic regulations.
- Make sure you are always visible to other road users. As a mobility scooter user, you are particularly vulnerable in traffic.
- Avoid quiet routes. Make sure help can be called in an emergency.

You may not use the mobility scooter if;

- you are physically unable to operate the vehicle safely due to, but not limited to, reduced muscle strength, reduced reaction time, or a lack of one or both arms
- you are unable to sit.

Assembly and adjustment

(see also the chapter 'Settings and adjustments')

Warning

- Before using the scooter, check that all components are properly assembled and securely fastened.

NB

- The dealer can use the microprocessor to set the various driving functions for you.

Intended user

The following levers may only be operated by the user:

- Throttle
- Brake lever
- freewheel lever
- lever for adjusting the seat depth
- lever for adjusting the backrest

The following levers may only be operated by the user and the assistant:

- lever for turning the chair
- hydraulic lever for the steering arm

Sitting in the scooter

Warning

- Avoid pinching your fingers and limbs. Always take extra care:
 - When folding away the armrests.
 - If you do not have full control over your arms or legs.
 - When you operate the freewheel lever

Driving the scooter

(see also the chapter 'Use')

**Warning**

- Drive carefully on roads that are slippery due to rain, icy conditions, or snow.
- Be careful not to get loose clothing caught on the wheels. This can cause the wheels to lock and the scooter to tip over.
- Prevent yourself from tipping over:
 - If you have adjusted the backrest to a large angle.
 - When driving up a slope.
 - When the sliding system is in the rear position.
- Do not carry anything on the handlebars.
- Do not allow anyone to stand on the back of the scooter.
- The scooter may suddenly come to a stop during use.

**NB:**

- In **the Netherlands** it is mandatory to insure your mobility scooter with at least third-party liability*.

* Third-party liability insurance (WA) is a liability insurance that is legally required by the government. Damage you cause to others with your mobility scooter, for which you are legally liable, is covered. This does not apply to damage caused to your mobility scooter through your own fault. An insured vehicle can be identified by an insurance plate on the rear.

Slopes and other obstacles**Warning**

- Avoid slopes with an incline of more than 10°. Driving on slopes (forward or backward) is always dangerous.
- Never freewheel down a slope. (See page 27: freewheel lever)
- Drive carefully along (pavement) edges. An uneven surface increases the risk of tipping.
- Avoid prolonged driving on slopes of less than

Maximum slope angle

straight from the front

Approaching a sidewalk

**NB**

- Your center of gravity shifts when you drive up or down a slope. Be careful.
- On a steep slope there is a greater risk of tipping over, so be careful when turning around.
- Always drive forward and in a straight line when driving on or off a sidewalk. Accelerate carefully. It's best to drive onto the sidewalk via a dip in the curb.
- The maximum height you can drive over is 100 mm.
- A soft, slippery or uneven surface can create a dangerous situation.

Transporting the scooter without a driver

(see also chapter 'Transport and storage')

**Life-threatening:**

- When transporting your mobility scooter in a motor vehicle, always ensure you are seated in a seat that cannot move. A mobility scooter should never be considered a seat in a motor vehicle, even if it is secured with a restraint system such as those used in a taxi.

**NB**

- The scooter must be fully assembled before use.

**Miscellaneous****Warning**

- Do not place the scooter in direct sunlight; the upholstery may become too hot and cause burns or skin irritation.

**Be careful**

- Prevent the scooter from coming into contact with seawater; seawater is corrosive and will damage the scooter.
- Prevent the scooter from coming into contact with sand, as this can penetrate the scooter's moving parts, causing them to wear out more quickly.

**NB**

A number of technical precautions have been taken to further ensure your safety:

- The immobilizer is activated when a battery charger is connected.
- When cornering, the speed is automatically limited.
- Freewheel is limited.
- A warning will appear on the display if the battery voltage is too low.
- Information appears on the display to help resolve possible faults.
- In the event of an electronic malfunction, the scooter will make a controlled stop.
- The scooter is equipped with an ergonomically shaped handlebar to prevent injuries to fingers and hands.
- To prevent body parts from becoming trapped, components where this could happen are shielded.
- Your dealer can install a seat belt (anterior hip support) on the armrest mounting points attached to the seatback. This will secure your body in the correct position.

Approval for EMC

The scooter is EMC (electromagnetic compatibility) approved. Tests show that a scooter is not sensitive to or harmful to mobile communication equipment and electronically controlled devices. Nevertheless, exercise caution when using the scooter simultaneously with mobile or electronic equipment.

**Warning**

- The driving characteristics of the scooter can be influenced by electromagnetic fields (from mobile telephones, power generators, high-voltage sources, etc.).
- The mobility scooter may disrupt the operation of devices in the immediate vicinity that emit electromagnetic fields (a shop's alarm system, automatic doors, etc.).

Introduction

The Solo RX is an electrically powered mobility scooter, suitable for young and old people who have difficulty walking long distances.

You can fully customize the Solo RX scooter to your liking. This guarantees a relaxed riding position. The scooter is operated with both hands. It's easy to operate, and you don't need to be strong to ride it. The Solo RX is available as a four-wheel scooter.

This manual contains the scooter's default settings. A list is included at the back of this manual where you can specify your specific settings.

Intended use

The Solo RX is designed as a mobility solution for people with limited mobility. It provides a seat support.

The user operates and controls the product from a seated position. The Solo RX is a product that can be configured based on the accompanying price list.

* Depending on the set speed, outside temperature, user weight, tire pressure, etc.

** Less is possible,

Technical specifications Solo 3 and 4

- Intended use: Outdoors, between -25° C and + 50° C, on paved paths and roads
- Dimensions (L x W x H):
 - Solo 3: 140 x 64 x 120 cm, (reduced: 140 x 64 x 90 cm)
 - Solo 4: 140 x 64 x 120 cm, (reduced: 140 x 64 x 80 cm)
- Height without chair 850 mm
- Weight without batteries:
 - Solo 3: 96 kg
 - Solo 4: 105 kg
- Weight with batteries (60 Ah): 96 kg
- Weight of heaviest removable part (seat): 20 Kg
- Maximum static stability downhill: 15°
- Maximum static stability uphill: 15°
- Maximum static stability lateral, uphill: 15°
- Maximum dynamic stability uphill: 10°
- Nominal slope angle: 10°
- Maximum theoretical driving distance according to ISO 7176-4 with fully charged batteries*: 100 km
- Maximum obstacle height 100 mm
- Maximum fall height obstacle (curb edge): 100 mm
- Ground clearance: 90 mm
- Max. (programmable) forward speed: 19 km/h**
- Minimum braking distance at maximum speed (programmable): 3600 mm
- Seating angle: 3°
- Seat depth: 440 mm
- Seat width: 510 mm
- Seat height above footrest (min. - max.): 500 - 580 mm
- Footrest height above ground: 80 mm
- Backrest angle: 90 - 160°
- Backrest height: 500 mm
- Armrest height (to seat): 220 mm
- Armrest length (to backrest): 300 mm
- Turning circle (radius)
 - Solo 3: 1400 mm
 - Solo 4: 1670 mm
- Turning space: 2000 mm
- Test dummy weight: 160 kg
- Maximum user weight: 160 kg

- Required operating force: about 1-30N
 - Battery capacity (optional): 60 Ah / 73 Ah / 80 Ah / 97 Ah / 160 Ah
 - Max. load: 8 amps
 - Engine, transaxle system 750 watt 24 V
 - Charger type: 8 Amp 24 V
 - Maximum noise level 75 Db
 - Magnetic brake: 10 Nm 24 V
 - Dismountable for transport: limited (see also page 32)
 - Wheel size: 3.00 - 8 (tubeless tire)
 - Complies with: EN 12184 class C
- The scooter has two 12-volt dry (gel) batteries. The following types are suitable:
- 60 Ah (198 x 166 x 171 mm) (14,2 kg)
 - 73 Ah (260 x 170 x 178 mm) (17,7 kg) (*Gel batterie*)
 - 97 Ah (260 x 170 x 210 mm) (23,6 kg) (*Gel batterie*)

Service and technical support

Consult your dealer if you have any questions about specific settings, maintenance, or repairs. Always have the serial number at hand. You can find this on the CE sticker on the scooter. (see also page 16).

Technical lifespan

Life & Mobility assigns this product a technical lifespan of 6 years. This means the product must be repairable, safe, and usable during this period.

The above declaration only applies when the user uses the product in the capacity for which it was intended by Life & Mobility and is CE marked as such, including any original components.

Product description


Inspection upon delivery

Before using your mobility scooter, please check that the delivery matches the product description below. If you notice any transport damage upon delivery, please report it to your dealer. Life & Mobility is confident that the mobility scooter will meet your expectations and hopes you will enjoy using it.


Product description

Mobility scooter

- A. Steering wheel with control panel
- B. Basket
- C. Battery charging contact
- D. Headlight
- E. Direction indicators (orange)
- F. Front wheel
- G. Lid CE sticker
- H. Foot platform
- I. Battery cover
- J. Freewheel lever
- K. Rear wheel
- L. Brake lights / Tail lights (red)
- M. Lumbar pump
- N. Lumbar release button
- O. Folding armrests
- P. Backrest
- Q. Ignition key
- R. Sliding carriage

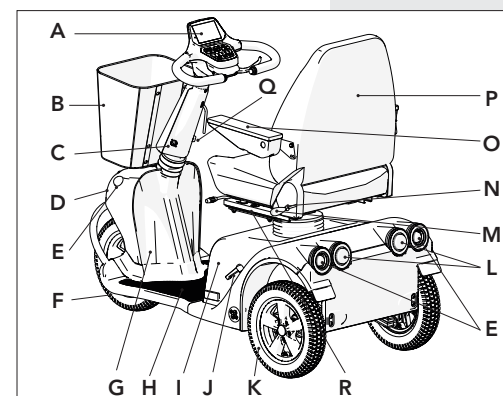Product description Solo 3

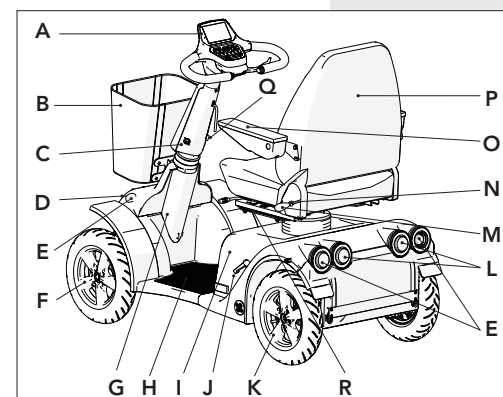
Product description Solo 4

Chair (sliding sled under the chair)


- S. Lever slide carriage
- T. Armrest adjustment wheel
- U. Back angle lever
- V. Chair lever

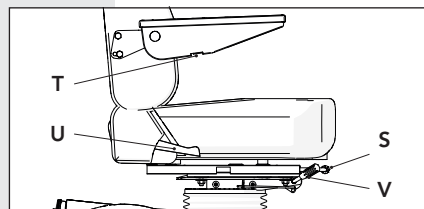Product description chair

Steering wheel


- W. Handlebar angle adjustment lever
- X. Allen screw handlebar height

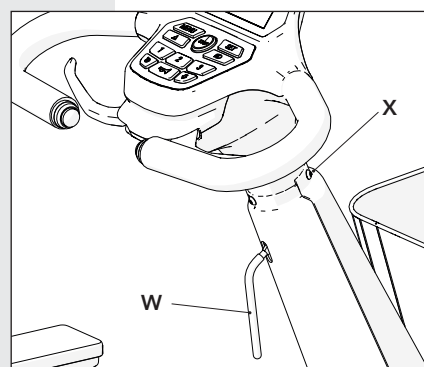Product description steering wheel

Control panel

- a. Emergency stop
- b. Hazard lights
- c. Battery indicator
- d. Speed level 1, 2, 3
- e. Throttle *
- f. Brake lever*
- g. Left/right turn signal*
- h. Cruise control
- i. Backward
- j. Horn
- k. Turn signal right/left*
- l. SET key
- m. Menu key
- n. Lighting
- o. Display

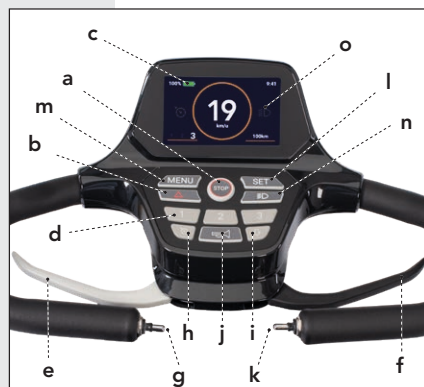Product description control panel

* Depending on how the scooter is programmed.

Suspension

You have a personal choice for optimal suspension comfort:

Active: Firmer suspension and a firm seat that allow for a sporty, active driving style.

Comfort: An optimal mix of shock absorption and driving behaviour.

Xtra: Extra soft suspension and soft seat, for (Only maximum shock absorption (e.g. due to pain caused by muscular disease, rheumatism or back problems).

Facilities

- Basket
- Lighting
- Horn
- Torsion suspension on front and rear wheels
- Semi-independent suspension on the rear axle
- Comfortable chair with pre-formed seat and backrest, with adjustable lumbar
- Simple and easy operation of the steering wheel and steering levers
- Speedometer, trip distance and total distance in km
- Clock and date
- Battery indication and low battery warning

Not pictured

Various accessories are available as options , including:

- Padded armrests
- Headrest
- Hip belt
- Additional foot throttle control
- Rearview mirrors
- Stick holders
- Rollator holder
- Extra shopping basket
- Lined lap mat, rain and wind resistant
- Rainproof cover

- Windshield
- 4 fixation eyes (Solo 3) / 2 fixation eyes (Solo 4)

For more accessories see www.life-mobility.com.

Indications on the mobility scooter

There are several stickers on the scooter. Below is an overview of their locations.


CE sticker

The Solo complies with the essential requirements of Regulation (EU) 2017/745 for Medical Devices and is provided with a CE mark.


The sticker for the maximum load capacity of the front and rear axle and the serial number are located under the cover (G) at the bottom of your steering column.

CE sticker

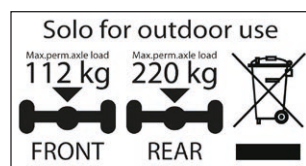
Axle load sticker

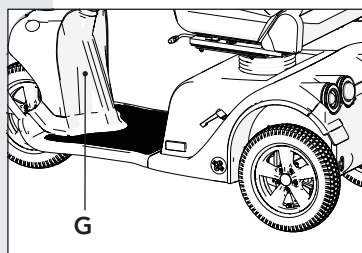A CE label is affixed to the steering column, indicating the maximum user weight and noting that it is a medical device [MD]. You are also reminded to read the manual.

The mobility scooter is not intended to be used as a seat in a motor vehicle.

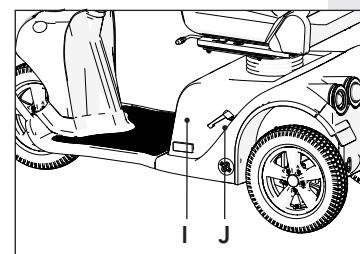Freewheel sticker (see also page 27)

The sticker above the freewheel lever (J) indicates the "push" position. In this position, the scooter can be moved freely.

The sticker under the freewheel lever indicates that the scooter is ready to drive electrically.

Battery connection diagram

The battery connection diagram is located on the inside of the battery cover (I). This also includes the fuse and engine control unit.

Push position

Electric driving mode


Settings and adjustments

You can configure the scooter entirely to your liking. We recommend consulting your dealer for the correct settings.

Adjusting the chair

Seat height

1. Move the lever (T) up (it is located on the left or right side of the seat depending on the selected setting).
When you sit in the chair, it lowers. When unloaded, it raises.
2. Release the lever (T).
The chair is locked in its current position.

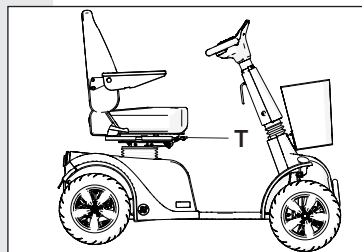
Change seat height

Seat depth

NB

- **YOU** Your center of gravity shifts when you move the chair using the slide. Be aware of this, especially when driving uphill.

1. Operate the lever (Q2)*.
Please note: do not confuse this with the seat height lever (T).
2. Slide the seat forward/backward.
3. Release the lever and the slide will click into place.

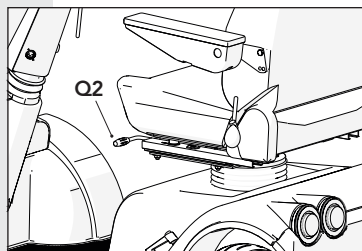Change seat depth


Back angle

NB

- Only adjust the backrest when you are sitting in the chair.
- The largest back angle (180°) is only intended for transporting the chair.

1. Take a seat in the chair
2. Pull the lever (S) up.

* This is located on the left or right, depending on how the scooter is adjusted.

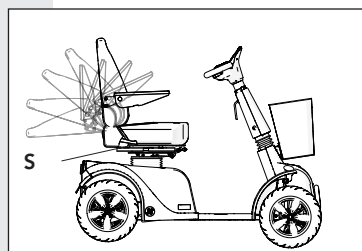Change back angle


3. Use your upper body to push the backrest back to the desired position, or:
4. Bend your upper body forward to make the backrest move forward automatically.
5. Release the lever (S).
The backrest is locked in its current position.

Armrests

NB

- The position of the armrests changes when the backrest is adjusted.


1. Rotate the armrest up/down to the desired angle.
2. Use the adjustment wheel (R) to fix the angle of the armrest.
3. Loosen the adjustment wheel (R) to move the armrest downwards.

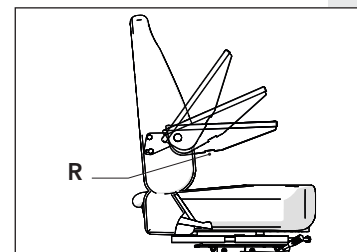Change armrest angle

Lumbar support backrest

1. Squeeze the lumbar pump (M) a few times.
The curvature of the backrest increases.
2. Check whether the new position has the desired result.
3. Press the lumbar release button (N) to reduce the curvature.

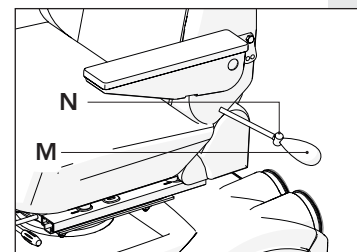Adjusting the lumbar backrest

Adjust steering wheel

Carefully


- Do not pull the handlebars up too hard when adjusting; the maximum adjustability is 5 cm.

NB

- Have the handlebar height properly adjusted by your dealer once in advance.
- Do not adjust the steering column angle while driving. The steering wheel can move forward automatically thanks to a gas spring.

Handlebar height

1. Loosen the Allen screw (X) with an Allen key.
2. Adjust the height of the handlebars.
3. Retighten the screw firmly.

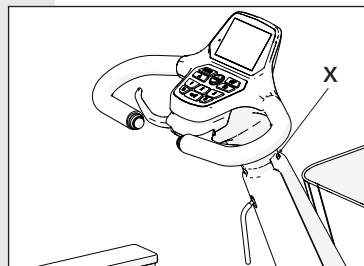
Adjust handlebar height


Steering column angle

1. Take a seat in the chair.
2. Pull lever (W) towards you.
3. Pull the steering column towards you or push it away from you until the desired position is reached.
4. Release the lever.


The steering column is now locked.

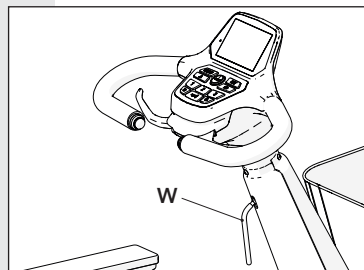Adjust steering angle

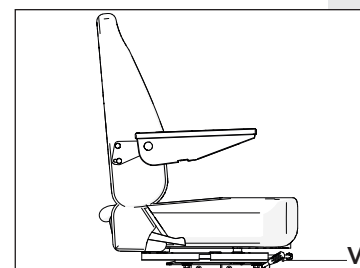Chair lever

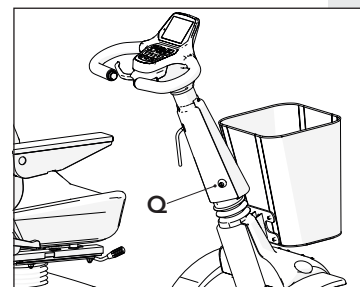Ignition key

Pre-check

Before driving away, check the following:

- the operation of the lighting;
- the battery voltage in the display;
- the tire pressure (see chapter 'Maintenance');
- whether the freewheel lever is pointing downwards (electric driving mode, see page 27);
- or the battery charger has been removed.
- Additional load may only be transported in basket B (maximum load: 3 kg)
- If you notice any unusual behaviour in your mobility aid (noise, odour, behaviour), please contact your dealer immediately.

Getting in and out

NB

- Make sure the seat is properly locked before driving.

Getting on / Getting off

1. Push the seat lever* (V) down.
2. Turn the chair a quarter turn to the left or right.
3. Get in or out.
4. Turn the chair back until you hear a clicking sound.

The chair is now locked.

Enable

1. Insert the key into the ignition (Q).
2. Turn the key a quarter turn clockwise.
3. Wait until '0 km/h' appears in the display.

The scooter is now ready to drive.

Maximum speed

The maximum speed can be set individually. Ask your dealer for advice.

** This is located on the left or right, depending on how the scooter is adjusted.*

**Warning**

- Control your speed before entering a corner. The cornering speed limiter will only adjust your speed if the corner is sharp. The scooter may tip over.
- Despite the scooter's small turning radius, turning into elevators or doors can be difficult. Always use the correct low speed setting for short maneuvers.

**NB**


- Control your speed with the throttle and brake levers. The speed setting determines how fast you're going when you fully depress the throttle.
- When reversing, check that the road behind you is clear.
- The scooter moves backwards more slowly than forwards.
- An audible signal will sound when reversing.
- Upon request, your dealer can change the way the throttle is operated.

The maximum speed is 19 km/h. Your dealer can adjust the speed setting.

1. Press button 1, 2 or 3 to select the speed setting.

Choose a suitable speed setting for each situation.

- Indoors, in shops and on the sidewalk: speed setting 1
- In busy streets: speed setting 2
- On quiet roads: speed setting 3

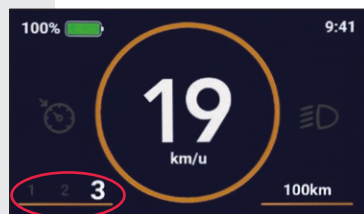Speed setting 1, 2 or 3

Driving forward and backward

1. Squeeze the right lever (f)* to drive forward.
2. To drive backward, push the "back key" button (i) and squeeze the right lever.

Switching between forward and reverse is only possible when the scooter has come to a complete stop.

Braking and reducing speed:**Remark:**

The lever on the left side is the brake lever.

1. Braking can be done by fully squeezing the brake lever or by releasing the throttle. Using the brake lever will bring the scooter to a stop more quickly.
2. You can reduce speed by partially squeezing the brake lever, which will cause the scooter to slow down. This can also be done by squeezing the throttle less deeply.

Cruise control

1. Drive at the desired speed.
2. Hold the throttle (e)* and press the cruise control button (h). You'll hear a beep. "CRUISE" will appear on the display.
3. Now release the throttle.
The cruise control is engaged.
4. To deactivate cruise control, operate the throttle (e)*, brake lever (f), or button (h). Cruise control automatically deactivates in sharp turns.

Lighting

NB

- After a sharp bend, the turn signals automatically switch off.
- The one-handed operation can be set up for you by your dealer.
- Your dealer can turn off the standard hazard warning lights chime for you.

1. Press key (n) to switch the lights on or off.
2. Use lever or push button (g) or (k) to switch the turn signals on or off.
3. Press the button (b) to switch the hazard lights on or off.

Batteries

If your scooter is equipped with lead-acid batteries:

Carefully

- New batteries need to be broken in. Batteries only reach their full capacity after several uses and charges. Keep this in mind during your first few rides.
- Charge the batteries daily if you use your scooter frequently, and every few days if you don't use them often. This will extend their lifespan.
- Always fully charge the batteries.
- Prevent the batteries from becoming completely discharged. This will cause the scooter to stall and negatively impact the battery's lifespan.
- Charge fully discharged batteries as soon as possible.
- Do not keep your batteries in a discharged state for extended periods of time.
- Never place metal parts on the batteries: this can cause a short circuit and damage the batteries.

NB

- At low temperatures, batteries have a lower capacity. This also reduces the range.

- Put the scooter on the charger every night.
- You cannot overcharge the included batteries. Thanks to voltage limitation, the charger automatically switches off when the batteries are full.
- Use only battery chargers from Life & Mobility.
- If you use batteries other than those supplied by Life & Mobility, ensure they are properly secured to prevent displacement. Failure to do so can lead to a short circuit.

If your scooter is equipped with lithium batteries:

Carefully

- Charge the batteries daily if you use your scooter frequently. This way, you'll have full range again the next day.
- Prevent the batteries from becoming completely discharged. The scooter will then stop without warning, and the display will also turn off.
- Charge batteries the same day if the charge is less than 30%.
- Never place metal parts on the batteries: this can cause a short circuit and damage the batteries.
- Lithium batteries have a built-in Battery Management System (BMS) that protects the battery. If the BMS detects improper battery usage, such as overcharging or excessive discharge, the battery will shut down and the scooter will no longer run. The display will also remain dark. The battery can be reactivated by connecting it to the charger. Depending on the cause, reactivation can take up to an hour.

NB

- At low temperatures, batteries have a lower capacity. This also reduces the range.
- You cannot overcharge the included batteries. Thanks to the voltage limiter, charging stops when the batteries are full.


- Use only Life & Mobility lithium battery chargers. Battery chargers for lead-acid batteries are not suitable.
- If you use batteries other than those supplied by Life & Mobility, ensure they are properly secured to prevent displacement. Failure to do so can lead to a short circuit.

Reading the battery status

The battery-shaped indicator on the screen shows the remaining battery capacity:


Battery status	Indication
100% Battery icon fully filled green	The batteries are full.
20% Battery icon partial fill is red	You can still drive. 20% of the total range (actual range depends on the battery type and other conditions). It's best to charge now.

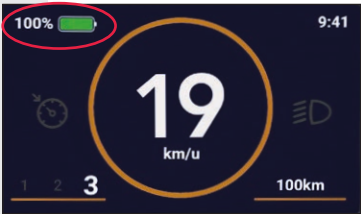Battery status indicator

Charging lead-acid batteries

1. Turn off the scooter with the ignition key.
2. Insert the battery charger plug into the battery charging socket (C).
3. Plug the power cord into a wall outlet.
4. Always complete the charging cycle. This can take up to 24 hours.

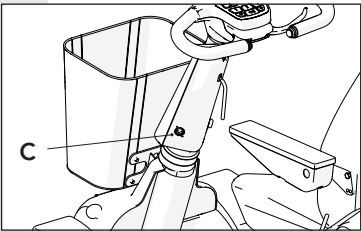Battery charging contact

Charging lithium batteries

1. Turn off the scooter with the ignition key.
2. Insert the battery charger plug into the battery charging socket (C).
3. Plug the power cord into a wall outlet.

NB


- The power plug must have been out of the wall socket for at least 1 minute to start a new charging cycle
4. It's best to complete the charging cycle. This can take up to 16 hours.
 5. Lithium batteries can also be partially recharged, for example, while traveling. In that case, make sure they are fully charged again within a week.

Freewheel lever

Warning

- Never use this function on a slope or while driving. The magnetic brake will then be inoperative. After each use of the freewheel function, turn the freewheel lever down to re-engage the brake.

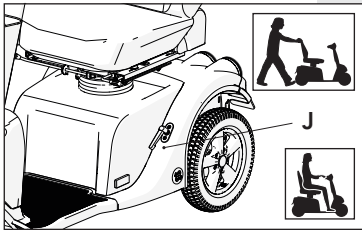Freewheel lever

1. To move the scooter manually, turn the freewheel lever (J) upwards. A beep sounds every 2 seconds to alert you that the freewheel is engaged.
2. Turn the freewheel lever back down for electric driving mode. Get used to this after each use.

NB

- If you accidentally roll down a slope in freewheel mode, the scooter will automatically slow down to a low speed. You can shorten the braking distance by pushing the freewheel lever down.

User menu

Setting the trip distance to 0

1. Ensure that the trip distance in kilometres (a) is visible, see 'Start screen'.
2. Press and hold the SET button (b) for 3 seconds. A beep will sound. The trip distance has been set to 0.

Setting the trip distance to 0

Setting the language

1. Press and hold the MENU button for 3 seconds.
(a)
2. Press the SET button (b) to enter the language change menu.
3. Press the MENU button (a) to select the next language in the list. Repeat this until the desired language is selected. Select 'Cancel' to not set a different language.
4. Press the SET button (b) to confirm.
5. Press the MENU button (a) for 3 seconds to exit the menu.

*Setting the language**Setting the time and date*

1. Press and hold the MENU button for 3 seconds.
(a)
2. Press the MENU button (a) to select 'Time and date'.
3. Press the SET button (b) to enter the Time and date change menu. In the change menu, the "hours" can be adjusted first.
4. Press the MENU button (a) to increase the hours by 1 hour. Repeat this until the desired number is displayed.
5. Press the SET button (b) to confirm the hours. You can now adjust the minutes.
6. Press the MENU button (a) to increase the minutes by 1 minute. Repeat this until the desired number is displayed and press the SET button (b).

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The date can now be adjusted

7. Use the SET button (b) to set the correct day, month and year, as well as the hours and minutes.
8. Press the SET button (b) to confirm and return to the user settings.
9. Press and hold the MENU button (a) for 3 seconds to return to the Home screen.

Cleaning*Carefully*

- Do not use bleach or solvents.
- Do not use a pressure washer to clean the scooter.

**NB**

- Clean the upholstery and frame regularly.

Upholstery

1. Remove light contamination with a solution of mild detergent and warm water.
2. Once the dirt has set into the upholstery, gently scrub the upholstery with a soft brush.

Frame and hoods

1. Clean the housing with a non-abrasive cleaning agent.
2. Apply a non-abrasive car wax to the entire frame and only to the painted fairings.

Batteries (see also chapter 'Use')**If your scooter is equipped with lead-acid batteries:****NB**

- Fully charge the batteries at least once a month.
- The battery connection method is illustrated in a diagram inside the battery cover. Note that this diagram is for lead-acid batteries.
- Be careful with damaged batteries. Leave the (environmentally friendly) replacement of batteries to your dealer.

*Carefully*

- Batteries can be dangerous. Have your batteries installed and maintained only by a qualified scooter mechanic. They have the necessary qualifications and tools to perform installation and maintenance safely and correctly.
- The scooter uses two 12-volt batteries with a 24-volt series connection. Never use simultaneously:

- batteries from different manufacturers;
- batteries with different technologies;
- Batteries with different capacities;
- batteries with different date codes.

**NB**

- Always fully charge batteries before first use.

If your scooter is equipped with lithium batteries:**NB**

- Fully charge the batteries at least once a month.
- The connection method of the batteries is explained on a diagram in the battery cover.
- Note that this is the diagram for lithium batteries.
- If the batteries haven't been charged for a long time, the battery indicator may fluctuate significantly. You then run the risk of the scooter stopping unexpectedly.
- Be careful with damaged lithium batteries. They can pose a fire hazard.
- Leave the (environmentally friendly) disposal of these batteries to your dealer.

**Carefully**

- Batteries can be dangerous. Have your batteries installed and maintained only by a qualified scooter mechanic. They have the necessary qualifications and tools to perform installation and maintenance safely and correctly.
- The scooter uses two 24-volt lithium batteries connected in parallel. Never use simultaneously:
 - Batteries with different voltages. If there is a voltage difference, fully charge both batteries separately before connecting them.
 - batteries from different manufacturers;
 - batteries with different technologies;
 - batteries with different date codes.

**NB**

- Always fully charge batteries before first use.

Tyres**Warning**

- For safety reasons, make sure your scooter's tire pressure is correct. Under-inflated tires will worsen handling and make steering harder.

**NB**

- The scooter's tires must have sufficient pressure. Check this regularly.
- Under-inflated tires will reduce your range.
- Too high tire pressure makes riding uncomfortable.
- Check the tires regularly for wear and (drying) cracks.

Service and technical support

Take your scooter to your dealer at least once a year for a service and inspection. Preventive maintenance involves checking the proper functioning of all moving parts, as well as the electronics and the condition of the batteries.

Repairs should be carried out exclusively by your dealer or by Life & Mobility, using only original parts supplied by Life & Mobility.

Maintenance matrix

The matrix below indicates the components that should be checked regularly. If you have any doubts or if a defect has occurred, you can contact your dealer to have it resolved.

Transport and storage

Element	Action	Daily	Weekly	Monthly	Per quarter	Per year
Battery charger	Check for operation	X				
Tyres	Check tire pressure (2.5 to 3 bar)		X			
Lighting	Check operation	X				
Housing	Clean with a damp cloth		X			
Wiring	Check for breakage		X			
Moving parts	Check for operation and play			X		
Batteries	Full charge cycle			X		
Anti-tip wheel (if available)	Check for operation				X	
Parking brake	Check for operation				X	
Tyres	Check for wear				X	
Mobility scooter	Service (by your dealer)					X
Carbon brushes	Check condition, replace if necessary (by your dealer)					X
Potentiometer	Check the potentiometer for play in the zero position*					X
Frame	Check for rust, deformations and other damage, carry out repairs if necessary.					X

* By moving the throttle, you can easily feel any play in the neutral position. Any unusual instability of the throttle in the neutral position requires extra attention (manufacturer's advice) or immediate replacement of the potentiometer.

Extremely dangerous:

Life-threatening

- When transporting your mobility scooter in a motor vehicle, always ensure you are seated in a seat that cannot move. A mobility scooter should never be considered a seat, even if it is secured with a restraint system such as those used in taxis.
- Secure the scooter in the car by running straps through the four rims and attaching them to the mounting points on the car floor.
- A driverless mobility scooter can be transported by plane, ship, or car. Ask the carrier if any special safety precautions are required.

Carefully

- Never lift the scooter by the plastic parts of the housing.

NB


- The chair is heavy. Ask for assistance if necessary.

Flatten the backrest

- Stand next to the scooter.
- Hold the backrest firmly.
- Pull the lever (S) up.
- Push the backrest forward.
- Release the lever (S).

Removing the chair

- Turn the lever (T) downwards.
- Turn the chair slightly.
- Lift the seat off the scooter.

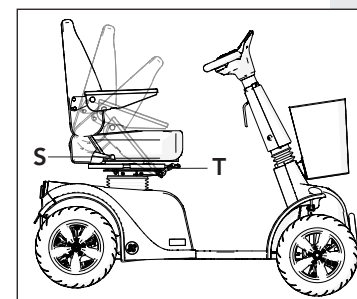Flatten the backrest

Laying the steering column flat

1. Stand next to the scooter.
2. Pull the angle adjustment lever (W).
3. Move the steering column as far as possible toward the seat.

Tip: When the seat is removed, the steering column can be laid further flat.



4. Release the lever again.

Transport

1. Turn the freewheel lever upwards.
2. Drive the scooter into the vehicle.
3. Turn off the scooter.
4. Make sure that the scooter is properly secured.

**Tip**

- Use 4 fixation eyes (Solo 3) / 2 fixation eyes (Solo 4).

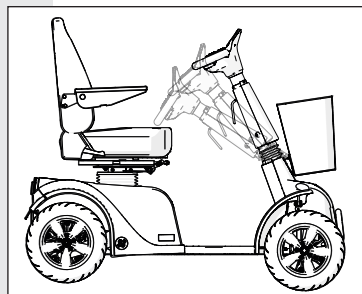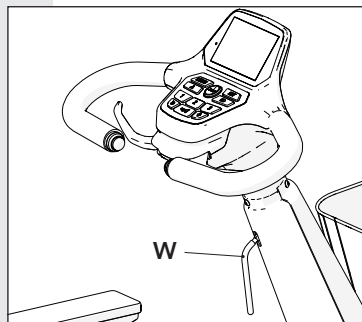Flatten the steering column

**Storage when not in use****NB**

- Always store batteries fully charged.
- Fully charged batteries can hold their charge for up to six months, but should always be charged before use.
- Avoid storing batteries at extremely high or extremely low temperatures.
- Sealed batteries do not require water refills. Consult a qualified scooter mechanic if you experience any performance issues.

Longer than a week

- Connect the battery charger at least once a week to maintain the battery's condition, or put the scooter away on the charger.

Very long lasting (e.g. during winter or due to illness)

- Have your dealer disconnect the battery wiring.
- After disconnecting the batteries, the battery indication may not be accurate until the scooter is fully charged again.

General**NB**

- Never try to fix a problem yourself. Contact your dealer.

List of possible error codes

1. The display reports a fault.
2. Press the menu key to read the fault text.
3. In many cases you can simply continue driving by turning the ignition off and on again once.
4. If the error persists on the display, please contact your dealer to have the fault resolved.
5. Report the fault to your dealer.

Warranty

Warranty conditions

Life & Mobility products are carefully assembled and meticulously inspected before leaving the factory. If a product does not meet your expectations, please contact the retailer from whom you purchased it.

Life & Mobility respects the laws and regulations applicable in the country of purchase to protect consumers.

The warranty does not cover defects caused by: failure to observe usage and maintenance instructions, use other than normal, wear and tear, carelessness, overloading, accidents involving third parties, the use of non-original parts, and defects caused by factors outside the product itself.

This warranty cancels all other warranties, whether stipulated by law or communicated orally, except those guaranteed in writing by Life & Mobility.

Guarantees are only valid within the EU.

Dispose of product and packaging

Dispose of product and packaging

All parts of your mobility scooter are recyclable. Your dealer is familiar with the NVMP collection points (ROS) for disposing of your mobility scooter.

From an environmental perspective, we recommend that you dispose of the packaging and the scooter separately as much as possible, so that the materials used can be recycled efficiently.

This product is suitable for reuse by another user. Your local dealer can prepare the scooter for another user. Before reuse, it must be cleaned and disinfected. Parts should be inspected for wear or breakage and replaced if necessary. Any necessary adjustments to size and functionality can be made according to the instructions in the manual.

This manual has been compiled with the utmost care, based on information known to Life & Mobility at the time of publication.

The information in this manual applies to the standard version of the product. Life & Mobility assumes no responsibility for any damage caused by changes or modifications to the product after the initial sale. Furthermore, Life & Mobility cannot be held liable for damage caused by work performed by third parties.

The information in this manual may not be reproduced and/or published in any way without the prior written permission of Life & Mobility.

The trade name and trademarks used by Life & Mobility may not be considered free under the Trade Name Act.

No rights can be derived from this text. Subject to change.

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© Life & Mobility – All rights reserved